

enziSurf[®] ROUTINE

multi enzymatic
surfaces
effect

• ECO-FRIENDLY • MEDICAL DEVICE •



enziSurf® ROUTINE

Multi-enzymatic decontaminant for cleaning of surfaces and floors in the healthcare sector (corridors, bedrooms, waiting rooms...) and care, medical, dental and veterinary practices).

- **MULTI-SURFACE** Cleans all surfaces without leaving a trace. Non-corrosive. Compatible with all materials.
- **FAST** Enzymes cut and disintegrate soil ensuring a faster cleaning with less effort.
- **DEEP CLEAN** The enzymes effectively and thoroughly remove visible and invisible organic dirt, even that which seems inaccessible. They also help to get rid of Coronavirus and harmful bacteria. enziSurf® fights odours.
- **PREVENTIVE** The anti-deposition effect of the enzymes prevents dirt from redepositing quickly.
- **NATURAL** enziSurf® offers a healthy alternative to conventional products, protecting the environment, staff and patients thanks to ingredients of natural origin (no phosphates, bleach or ammonia).


Instructions for use

Cleaning is an essential step in infection control.

1. Dilute the product in **hot water** as follows:
 - **2x30ml** in **5L of warm water** or
 - **120ml** in **10L of hot water**

Recommended operating temperature between 35-50°C


2. Apply the product to the surface to be decontaminated.
3. Scrub if necessary
4. Let dry.

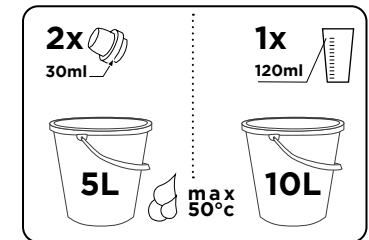

NOTE: make sure to regularly clean the surfaces in contact with patients. In addition to cleaning, it is recommended to disinfect in order to finalize the decontamination process. Nonetheless, the Superior Health Council recommends using disinfection sparingly in order to fight anti microbial resistance (AMR).

Applications

- Tables
- Working table
- Door handles
- Offices
- Counters
- Floors
- Corridors

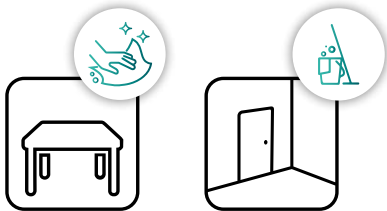

Product Characteristics

- Concentrated formulation combining raw materials of plant origin with enzymes for greater efficiency
- Non-aggressive for the user and the material in contact (compatible with all materials)
- Odourless
- Dye free
- Highly active enzyme complex
- Rate of biodegradability >97% in 28 days (OECD 302B)



Chemical Properties

Appearance: yellowish liquid

Precautions

Store in the original container, closed, between +4°C and +25°C

End of life: 24 months after production

Main Components

99% of the total ingredients are of natural origin.

Enzymatic complex:



protease
amylase
lipase

Delivery Units:



Ref **OL52632**
6 x 1L



Ref **OL52633**
1 x 5 L



Enzymes for patients' safety

Contact:

OneLife s.a.
Avenue Albert Einstein, 15
1348 Louvain-la-Neuve (Belgium)
+32 10 48 34 27 info@onelif-bf.com

onelif-bf.com